

Wax Eater

Ear Speculum/Otoscope Cleansing device


Any ear examined with an otoscope is going to result in a contaminated ear cone speculum. The **Wax Eater** otoscope speculum cleaner is designed to mechanically remove wax plus disinfect the otoscope speculum between each and every clinical use. This reduces the risks of inadvertently spreading pathogenic organisms i.e. yeast, mites or even *Pseudomonas* spp. to patients.



The specially designed brushes get into hard to reach areas of the otoscope speculum. The brushes are in a disinfectant solution of choice. A **Wax Eater** should be readily available in every exam room or anywhere where an otoscope is used. Weekly cleaning is recommended and we suggest attaching a piece of tape with last cleaning date to base of **Wax Eater**. Brushes should be replaced after showing signs of heavy wear.

- Compatible with most antimicrobial disinfectants
 - Rapidly cleans speculum
 - Works on most sizes of speculums
 - Limited one year warranty
 - 3 month on brush component
- Size: 5" (13.5cm) high
Diameter: 1¼"
Weight: 5oz

J1184 Wax Eater
J1184D1 Wax Eater replacement brushes
Brushes are compatible with other brands

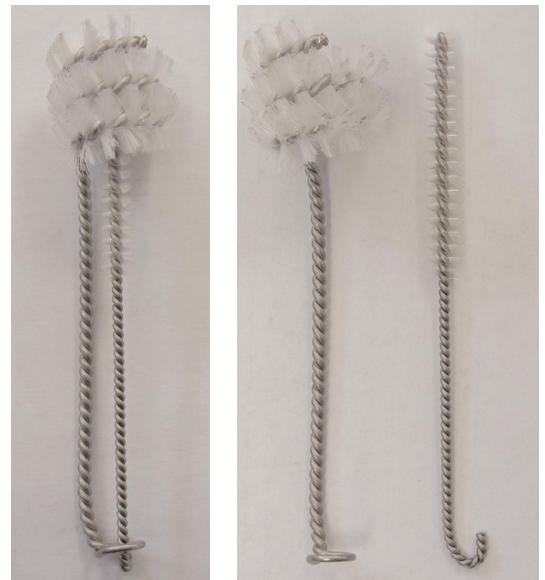

Jorgensen Laboratories, Inc.
Loveland, CO
(800) 525-5614 www.JorVet.com Info@JorVet.com