



Veterinary Energy-based Systems





## Veterinary Energy-based Systems

*Covidien is the world leader in designing and manufacturing advanced energy-based medical systems and accessories worldwide. Our product offerings include the complete line of Covidien state-of-the-art electrosurgical generators and an extensive array of electrosurgical accessories. Whether your surgical requirements call for a cutting edge vessel sealing system for permanent vessel occlusion or for the perfect low power office generator with monopolar/bipolar capability, look no further than Covidien Animal Health.*



### ForceTriad™ Energy Platform

- Valleylab™ mode for a unique combination of monopolar hemostasis and dissection while using a lower power setting
- Instant Response™ technology for precise energy delivery to effectively seal the vessel or tissue bundle
- Next generation LigaSure™ system utilizes TissueFect™ technology for sealing vessels 2x faster than previous technology
- Bipolar instrument receptacle for a traditional bipolar instrument
- Monopolar output equipped with Return Electrode Monitoring (REM™) for reduced risk of patient burns

## LigaSure™ Vessel Sealing System

- Instant Response™ technology for precise energy delivery to effectively seal the vessel or tissue bundle
- LigaSure™ tissue fusion technology for vessels and tissue bundles up to and including 7 mm
- Bipolar instrument receptacle for a traditional bipolar instrument




## Force FX™ Electrosurgical Generator

- Instant Response™ technology for delivering consistent power as tissue impedance changes
- Monopolar output equipped with Return Electrode Monitoring (REM™) for reduced risk of patient burns
- Bipolar instrument receptacle for a traditional bipolar instrument



### Tissue Response Technology vs Conventional Technology

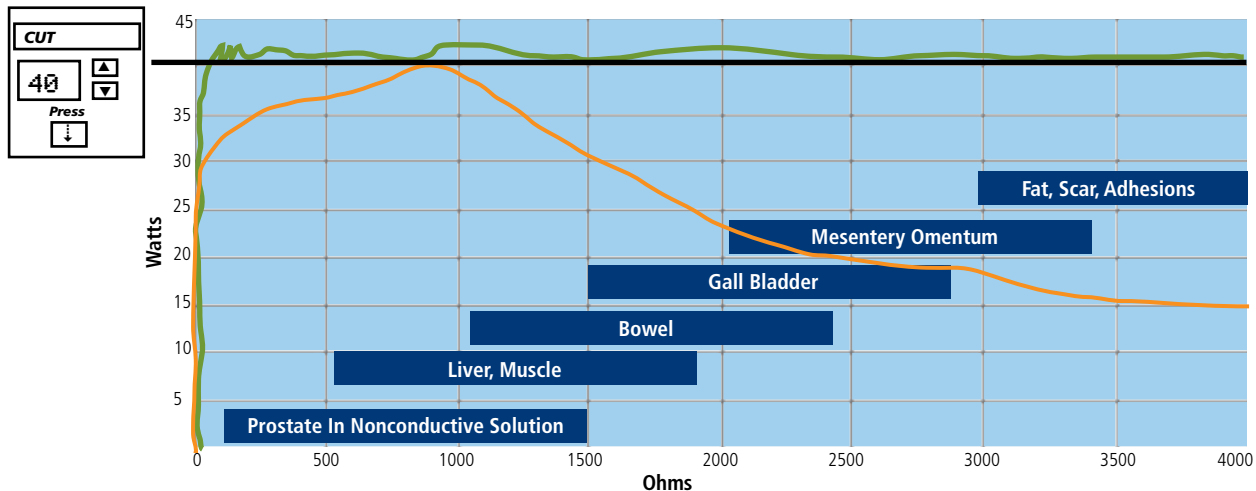

Tissue Response Technology™ senses tissue impedance at the electrode site and adjusts current and output voltage 200 times per second to maintain a consistent surgical effect. This is the principle behind Covidien Instant Response™ technology.

- Tissue Response Technology
- Conventional Technology



## SurgiStat™ II Electrosurgical Generator

- Bipolar instrument receptacle for a traditional bipolar instrument
- Monopolar output equipped with a return electrode contact quality monitoring system (RECCMS) for increased patient safety

### Accessories



Force Triverse™ Pencil



REM Polyhese™ II Patient Return Electrode



Ligasure™ Atlas



Buttonswitch Pencil



Bipolar Forceps

### Order Information

Code	Description	Code	Description
FORCEFXC	Force FX™ FX-C Generator	E751025	Infant Return Electrode
FORCETRIAD	ForceTriad™ Energy Platform	LS0300	LigaSure™ Single Footswitch
SURGI20	SurgiStat™ II Generator	LS1037	LigaSure™ Atlas Laparoscopic SLR/DIV
LIGASURE	LigaSure™ System Generator-Dom	LS1020	LigaSure™ Atlas 20cm
FT3000	Force Triverse™ Electrosurgical Device 3m (10Ft)	LS1200	LigaSure™ Precise
E6008	Standard Footswitch	LS1500	LigaSure™ V Sealer/Driver
E2250H	Footswitch Pencil	LS3090	LigaSure™ Max Instrument
E2350H	Rockerswitch Pencil	LS 3091	LigaSure™ Max Disp Elect/Cord
E2450H	Buttonswitch Pencil	LS3110	LigaSure™ TM XTD Instrument
E7507	REM Polyhese™ II Patient Return Electrode	LS3111	LigaSure™ TM XTD Disp Elec/Cord



For more information contact your local Covidien Animal Health sales representative or e-mail us at: [animalhealth@covidien.com](mailto:animalhealth@covidien.com)



COVIDIEN, COVIDIEN with Logo and ™ marked brands are trademarks of Covidien AG or an affiliate.  
©2008 Covidien AG or an affiliate. All rights reserved.

H6338 10M 09/08

15 HAMPSHIRE STREET  
MANSFIELD, MA  
02048

1-800-962-9888  
508-261-8000

[WWW.COVIDIEN.COM](http://WWW.COVIDIEN.COM)