

scil Vet abc Plus+ Hematology



About scil Vet abc Plus+

The scil Vet abc Plus+ Hematology Analyzer is the next level in hematology testing. From only 10 µL of blood, the scil Vet abc Plus+ performs a complete blood count including a 4-part differential of white blood cells. It provides you with essential information for fast diagnosis and treatment. The scil Vet abc Plus+ identifies all parameters of a complete blood cell count: erythrocytes, thrombocytes, leukocytes, hemoglobin and all erythrocyte-indices. A 4-part differential is available for the dog, cat, ferret and horse. Lymphocyte percentage/absolute number, Monocytes percentage/absolute number, Granulocytes percentage/absolute number and Eosinophil percentage/absolute number.

Features & Benefits

- \$2.95 cost per test (based on 4-5 test per day)
- 18 parameters in 60 seconds
- Small sample size - 10 µL
- Requires no sample preparation, no spinning
- Intuitive touch screen
- Automated QC +/- calibration
- Levy-Jennings charts on QC
- Compatible with most practice management software systems
- 4-part differential
- Stores up to 1,000 patient results
- 2-year warranty, 24/7 support

Data Management

Internal Data Storage: up to 1,000 results


Computer Interface: RS232, Ethernet

Sample Size: 10 µL EDTA blood

Test Time: 60 seconds

Throughput: 50 samples per hour

Report of Results: Display via touchscreen, printout, transfer to practice management programs

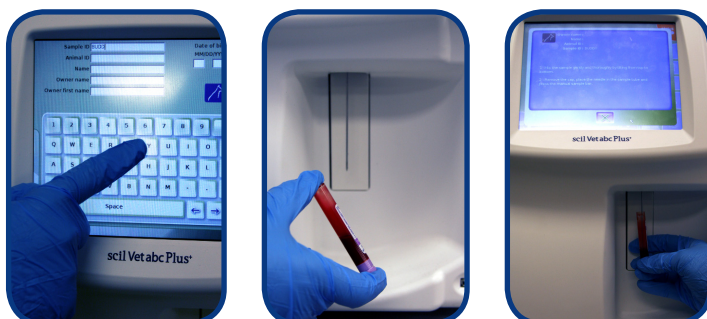

Parameters

WBC, RBC, HCT, HGB, PLT, MCV, MCH, MCHC, RDW, MPV, GRAN (# and %), MON (# and %), LYM (# and %), EOS (# and %)

Units: conventional or SI

Principle of Operation: Impedance technology for WBC, RBC, PLT. Spectrophotometry for hemoglobin

Reagents: Self-contained reagent pack system including diluent, lyse, cleaner and waste collection.

4-Part Differential: Granulocyte, Monocyte, Lymphocyte, Eosinophil

Animal Species

- Dog
- Cat
- Horse
- Rabbit
- Mouse
- Rat
- Sheep
- Pig
- Cattle
- Ferret

 HENRY SCHEIN®



scil animal care
a division of Henry Schein Animal Health
151 N. Greenleaf Street, Gurnee, IL 60031
P 1.877.724.5838 | www.scilvet.com